


Alere™ iCassette® Procedure Card



Test in 3 Easy Steps

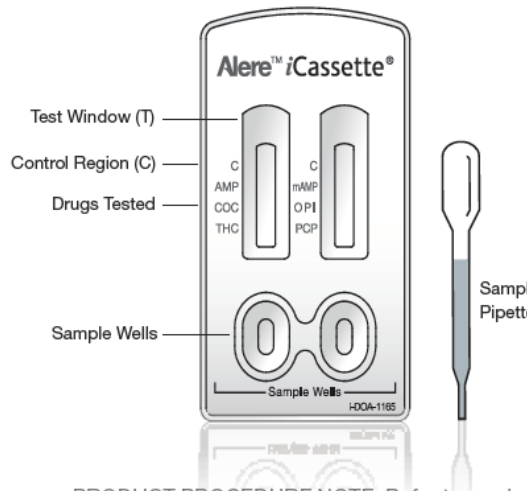

1
Remove the test from the sealed pouch.

2
Add 3 drops of urine to each specimen well.

3
Interpret results at 5 minutes.

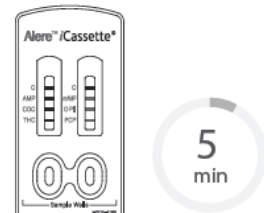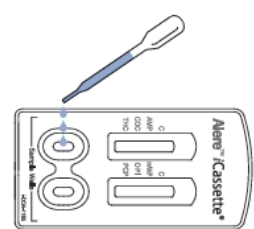

PRODUCT PROCEDURE NOTE: Refer to product insert for complete instructions, limitations, and warnings. For healthcare professionals including professionals at point of care sites. For in vitro diagnostic use only. This assay provides only a preliminary analytical test result. A more specific alternate chemical method must be used in order to obtain a confirmed analytical result. Gas chromatography-mass spectrometry (GC-MS) is the preferred confirmatory method. Clinical consideration and professional judgment should be applied to any drug of abuse test result, particularly when preliminary positive results are indicated.

Made by:

